



UNITY V3 AND UNITY VIA II/VR PROGRESSIVE LENSES								
Design	VSP Code	Compensated	CR-39	Polycarbonate	Trivex	1.60	1.67	1.74
Unity V3 Elite/Unity Via Elite II	N	•	•	•	•	•	•	•
Unity V3 Elite Fixed/Unity Via Elite II Fixed	N	•	•	•	•	•	•	•
Unity Via Elite VR	N + TA	•	•	•	•	•	•	•
Unity V3 Wrap	N	•	•	•	•	•	•	•
Unity Via Wrap II	O	•	•	•	•	•	•	•
Unity V3 Mobile	N	•	•	•	•	•	•	•
Unity Via Mobile II	O	•	•	•	•	•	•	•
Unity V3 Plus/Unity Via Plus II	O	•	•	•	•	•	•	•
Unity V3/Unity Via II	F		•	•	•	•	•	•
Lens Options			All SunSync Products/ Polarized	All SunSync Products/ Polarized	All SunSync Products/ Polarized	Polarized	All SunSync Products/ Polarized	

**Patient Benefit:** Clearer vision near, far, and everywhere in-between with a single pair of glasses.

**Visual Need:** Presbyopes requiring multifocal lenses.

UNITY V3 POWER RANGES V3 ELITE/V3 PLUS/V3 MOBILE/V3			
MATERIAL	SPHERE	CYL	ADD
Plastic 1.50	-13.25 to +6.50	To -6 [D]	+0.50 - +4.50
Poly 1.59	-15.00 to +9.00	To -8 [D]	+0.50 - +4.50
Trivex 1.53	-13.50 to +7.00	To -7 [D]	+0.50 - +4.50
Hi-Index 1.60	-15.00 to +9.00	To -8 [D]	+0.50 - +4.50
Hi-Index 1.67	-18.00 to +10.00	To -9 [D]	+0.50 - +4.50
Hi-Index 1.74	-19.00 to +12.00	To -10 [D]	+0.50 - +4.50

UNITY V3 WRAP POWER RANGES			
MATERIAL	SPHERE	CYL	ADD
Plastic 1.50	-7.00 to +6.00	to -4.0 [D]	+0.50 - +4.50
Poly 1.59	-6.50 to +6.00	to -5.0 [D]	+0.50 - +4.50
Trivex 1.53	-7.00 to +7.00	to -6.0 [D]	+0.50 - +4.50
Hi-Index 1.60	-7.00 to +7.00	to -6.0 [D]	+0.50 - +4.50
Hi-Index 1.67	-8.00 to +8.00	to -7.0 [D]	+0.50 - +4.50
Hi-Index 1.74	-9.00 to +9.00	to -8.0 [D]	+0.50 - +4.50

UNITY VIA OFFICEPRO LENSES						
Design	VSP Code	CR-39	Polycarbonate	Trivex	1.60	1.67
Unity Via OfficePro 5 FT.	IA	•	•	•	•	•
Unity Via OfficePro 10 FT.	IA	•	•	•	•	•
Lens Options		All SunSync Products/ Polarized	All SunSync Products/ Polarized	All SunSync Products/ Polarized	Polarized	All SunSync Products/ Polarized

**Patient Benefit:** Improved visual performance, as well as reduced eyestrain, fatigue, neck, back, and shoulder pain when working on a computer or a wide variety of workplace settings.

**Visual Need:** For progressive wearers who work on a computer for more than two hours per day, or perform tasks requiring extensive use of intermediate vision.

ALL UNITY VIA II DESIGNS POWER RANGES					
Material	Power Range	Power Range Unity Via II Wrap	Add Power Unity Via II Designs	Add Power Unity Via OfficePro	Cylinder*
CR-39	-7.00 to +6.00	-7.00 to +6.00	+0.50 to +4.50	+1.00 to +3.50	Up to -6
Polycarbonate	-10.00 to +6.50	-6.50 to +6.00	+0.50 to +4.50	+1.00 to +3.50	Up to -6
Trivex	-10.00 to +6.00	-5.00 to +4.00	+0.50 to +4.50	+1.00 to +3.50	Up to -6
1.60	-10.00 to +6.00	-6.50 to +5.25	+0.50 to +4.50	+1.00 to +3.50	Up to -6
1.67	-12.00 to +8.25	-6.00 to +5.00	+0.50 to +4.50	+1.00 to +3.50	Up to -6
1.74	-15.00 to +10.00	-	+0.50 to +4.50	+1.00 to +3.50	Up to -6

\*Prism in Rx's are dependent on the lab's capabilities. Unity V3 and UVII designs accommodate wide ranges of prism values. Please contact your lab to confirm availability.

UNITY RELIEVE							
Design	VSP Code	CR-39	Polycarbonate	Trivex	1.60	1.67	1.74
Unity Relieve 50	BA + TA	•	•	•	•	•	•
Unity Relieve 70	BA + TA	•	•	•	•	•	•
Lens Options		All SunSync Products/ Polarized	All SunSync Products/ Polarized	All SunSync Products/ Polarized	Polarized	All SunSync Products/ Polarized	

**Patient Benefit:** Single vision lenses that reduce the effects of single vision lenses that help reduce exposure to high-energy blue light. From device use.

**Visual Need:** For single vision patients who use digital devices for more than two hours per day.

### Reminders

UNITY SINGLE VISION							
Design	VSP Code	CR-39	Polycarbonate	Trivex	1.60	1.67	1.74
SVx, SVxtra, SVxtreme	BA	•	•	•	•	•	•
Lens Options		All SunSync Products/ Polarized	All SunSync Products/ Polarized	All SunSync Products/ Polarized	Polarized	All SunSync Products/ Polarized	

PREFERRED ENHANCEMENTS			
Light-Reactive Lenses	VSP Photochromic Category	VSP Code	Description
SunSync Elite (gray/brown)	Photochromics -Plastic	PR	<ul style="list-style-type: none"> <li>• Ultra-fast reaction times</li> <li>• UV protection and targeted blue light filtration</li> <li>• Excellent indoor clarity and outdoor comfort</li> </ul>
SunSync (gray/brown) SunSync Drive XT (gray/brown)	Photochromics -Plastic	PR	<ul style="list-style-type: none"> <li>• Fast reaction times</li> <li>• UV protection and helps filter high-energy blue light</li> <li>• Excellent indoor clarity and outdoor comfort</li> <li>• Retains color in the car (SunSync Drive XT)</li> </ul>
Anti-Reflective Coatings	VSP AR Coating Category	VSP Code	Description
TechShield® Elite (optional UV <sub>R</sub> )	D	QV (QV + BV)	<ul style="list-style-type: none"> <li>• Superior visual performance; enhanced cosmetic appearance</li> <li>• Improved scratch resistance, smudge resistance, cleanability</li> <li>• Super oleophobic; super hydrophobic</li> <li>• Optional backside UV protection</li> </ul>
TechShield Blue	D	QV (QV + BV)	<ul style="list-style-type: none"> <li>• Targeted blue light defense</li> <li>• Helps minimize exposure to high-energy blue light.</li> <li>• Superior visual performance; enhanced cosmetic appearance</li> <li>• Improved scratch resistance, smudge resistance, cleanability</li> <li>• Super oleophobic; super hydrophobic</li> <li>• Backside UV protection included</li> </ul>
TechShield Plus (optional UV <sub>R</sub> )	C	QT (QT + BV)	<ul style="list-style-type: none"> <li>• Superior visual performance; enhanced cosmetic appearance</li> <li>• Improved scratch resistance, smudge resistance, cleanability</li> <li>• Oleophobic; hydrophobic</li> <li>• Optional backside UV protection</li> </ul>
TechShield (optional UV <sub>R</sub> )	C	QT (QT + BV)	<ul style="list-style-type: none"> <li>• Superior visual performance; enhanced cosmetic appearance</li> <li>• Improved scratch resistance, smudge resistance, cleanability</li> <li>• Hydrophobic</li> <li>• Optional backside UV protection</li> </ul>
UV <sub>R</sub>		BV	<ul style="list-style-type: none"> <li>• Reduces back surface reflections</li> </ul>

### Reminders




- Pair all SunSync products with TechShield Elite (not TechShield Blue) for optimal light-reactive performance.
- Charge additional \$10 for backside UV on any TechShield AR Coating with UV<sub>R</sub>.

**Visual Need:** Ideal for prescriptions exceeding +3.00 or -3.00 and/or greater than 1.50 cylinder.

ALL UNITY SINGLE VISION DESIGNS—POWER RANGES		
Material	Power Range	Cylinder*
CR-39	-12.25 to +6.00	Up to -6
Polycarbonate	-14.00 to +6.00	Up to -6
Trivex	-14.00 to +6.00	Up to -6
1.60	-14.00 to +6.00	Up to -6
1.67	-16.00 to +8.00	Up to -6
1.74	-16.00 to +10.00	Up to -6

\*The cylinder range is an estimated range that can be produced properly, but in terms of calculus, the range is wider. Listed are basic availability guidelines. Upper limits may vary based on frame size, decentration, and blank size.

## UNITY DIGITAL PRODUCT PORTFOLIO

				
PRODUCT PORTFOLIO	DIGITAL FREEFORM PROGRESSIVE	DIGITAL FREEFORM SINGLE VISION	ANTI-REFLECTIVE COATINGS	LIGHT-REACTIVE LENSES
<b>BEST</b>	<b>Unity V3 Elite - Variable and Fixed</b>  <b>Unity V3 Wrap</b>  <b>Unity V3 Mobile</b> <i>VSP N</i>	<b>Unity Relieve</b>  50 or 70 Recommend pairing with TechShield Blue  <i>VSP BA + TA (\$10)</i>	<b>TechShield Blue</b> (for blue light reduction)  Backside UV included (add'l \$10 unless UV Covered)  <u><i>D + BV (mandatory)</i></u> <b>TechShield Elite</b> Optional Backside UV (\$10)  <i>D + *optional BV</i>	<b>Ultra Fast</b>  <b>SunSync Elite</b> Dark to clear in seconds. <i>VSP Category</i> Photochromics - Plastic  <i>VSP PR</i>
	<b>Unity Via Elite II - Variable and Fixed</b>  <i>VSP N</i>  <b>Unity Via Elite VR</b> <i>VSP N + TA (\$40)</i>			<b>Ultra Dark</b>  <b>SunSync Drive XT</b> Works behind the windshield. <i>VSP Category</i> Photochromics - Plastic  <i>VSP PR</i>
	<b>BETTER</b>			<b>Unity V3 Plus</b> <b>Unity Via Plus II</b> <b>Unity Via Wrap II</b>  <i>VSP O</i>
<b>GOOD</b>	<b>Unity V3</b> <b>Unity Via II</b>  <i>VSP F</i>	<b>Unity SVx</b>  <i>VSP BA</i>	<b>TechShield</b>  Optional Backside UV (\$10)  <i>VSP C</i>	
<b>LIFESTYLE</b>	<b>Unity Via OfficePro</b> (for computer/workplace/task)  5 ft. and 10 ft.  <i>VSP IA</i>  <b>Unity Via Mobile II</b> (for handheld digital devices)  <i>VSP O</i>	<b>Unity SVxtreme</b> (for wrap frames)  <i>VSP BA</i>	<b>MAXIMIZE PATIENT EXPERIENCE</b> <ul style="list-style-type: none"> <li>Discuss patient's visual needs for everyday, work, and play.</li> <li>Recommend solutions to address their needs.</li> <li>Explain how recommendations will provide them their best vision experience.</li> </ul>	

**Notes**

- Position of Wear measurements optional with O and N Category Unity Via Designs
- Simply take Monocular Distance PD and OC height for ALL Unity Designs
- \$10 copay for optional backside UV on TechShield Elite, TechShield Plus, TechShield
- SunSync products available in both Gray and Brown

**Guarantees**

- All Unity lenses backed by a satisfaction guarantee
- SunSync Light-Reactive Lenses one-year, 100% satisfaction guarantee (VSP members)
- TechShield AR Coatings two-year warranty

# UNITY QUICK REFERENCES

## Fitting Steps


1. Adjust the frame before measuring.
2. Measure monocular PDs, OCs, and/or fitting heights.
3. Fit from the center of the pupil to the deepest part of the lens.
4. Take Position of Wear (POW) measurements for a personalized visual solution (optional for category O and N lenses).

**Note:** Defaults are used when POW is not taken.

All Unity Lenses					
Design Name	Minimum Fitting Height	Personalized	Default Panto	Default Wrap	Default BVD
Unity SVx	N/A	No	-	-	-
Unity SVxtra	N/A	Yes	8	5	14
Unity SVxtreme	N/A	Yes	8	15	14
Unity V3 Plus/Unity Via Plus II	Variable length. Shortest = 12 mm	Yes	12	5	14
Unity V3 Mobile/Unity Via Mobile II	Variable length. Shortest = 12 mm	Yes	12	5	14
Unity V3 Wrap/Unity Via Wrap II	Variable length. Shortest = 12 mm	Yes	12	15	14
Unity V3 Elite/Unity Via Elite II	Variable length. Shortest = 12 mm	Yes	12	5	14
Unity Via OfficePro 5	Variable length. Shortest = 14 mm	Yes	12	5	14
Unity Via OfficePro 10	Variable length. Shortest = 14 mm	Yes	12	5	14
Unity Via Elite VR	Variable length. Shortest = 12 mm	Yes	12	5	14

## Unity Via OfficePro

- Minimum OC height of 14mm.
- Vertical distance between OC and upper inside edge of frame should be at least 10mm.
- Once you receive the patient's Unity Via OfficePro Lenses, remember to verify total near power in the reading area ONLY.

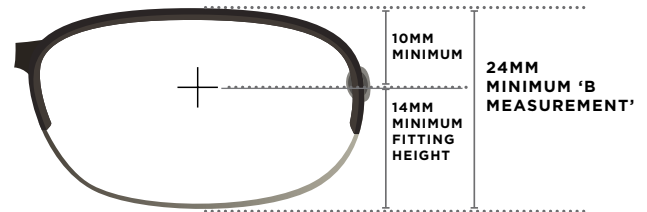

## All Unity Lens Designs

- Pre-adjust frame before measuring for optimal customization and performance.
- Ensure patient's eyes are centered in the lens (even more important if the patient's distance RX is greater than +/-4.00).
- Take monocular distance PD and OCs.
- Provide lab actual OC measurement as taken; DO NOT drop from measured OC.
- Ideal Position of Wear fit parameters:
  - Back vertex distance -14 mm by adjusting the nose pad
  - Pantoscopic tilt -12°
  - Positive wrap angle (bends toward the face) of around -5°
  - Default measurements are used if none provided

Learn more from your Ophthalmic Business Consultant, or visit [unitylenses.com](https://www.unitylenses.com).



## BEST PRACTICES:

### Maximize Patient Experience

- Discuss patient's visual needs for everyday, work, and play.
- Recommend solutions to address their needs.
- Explain how recommendations will provide them their best vision experience.

### Maximize Member Benefits

- Verify VSP member coverage prior to patient visit.
- Maximize VSP member benefits so patients receive the greatest value and the best vision experience possible (e.g., if progressives are covered, offer N category).
- Apply special member offers, multiple pair discounts, and/or rebates.

### Ordering Lenses

- Pair Unity Lenses with Preferred Enhancements.
  - Choose SunSync Light-Reactive Lenses for photochromic orders.
  - Choose TechShield AR Coatings for AR orders.
- Select Unity by name in the alphabetical drop-down menu when ordering.
- **Don't** select "Lab Choice" when ordering Unity products to ensure you receive the Unity lens you've recommended for your patient's best visual experience.

### Guarantees

For VSP members, contact VSP Insurance for a new authorization for Unity, SunSync, and TechShield guarantees.

Learn more from your Ophthalmic Business Consultant, or visit [unitylenses.com](https://www.unitylenses.com).

